

Evaluation of Two Methods for the Detection of Third Generation Cephalosporins Resistant Enterobacterales Directly From Positive Blood Cultures

Clarisse Durand, Agathe Boudet, Jean-Philippe Lavigne, Alix Pantel

▶ To cite this version:

Clarisse Durand, Agathe Boudet, Jean-Philippe Lavigne, Alix Pantel. Evaluation of Two Methods for the Detection of Third Generation Cephalosporins Resistant Enterobacterales Directly From Positive Blood Cultures. Frontiers in Cellular and Infection Microbiology, 2020, 10, pp.491. 10.3389/fcimb.2020.00491. hal-03334166

HAL Id: hal-03334166 https://hal.umontpellier.fr/hal-03334166v1

Submitted on 3 Sep 2021

HAL is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers. L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.



Distributed under a Creative Commons Attribution 4.0 International License



Evaluation of Two Methods for the Detection of Third Generation Cephalosporins Resistant Enterobacterales Directly From Positive Blood Cultures

Clarisse Durand¹, Agathe Boudet², Jean-Philippe Lavigne^{2*} and Alix Pantel²

¹ Service de Microbiologie et Hygiène Hospitalière, CHU Nîmes, Nîmes, France, ² VBMI, INSERM U1047, Université de Montpellier, Service de Microbiologie et Hygiène Hospitalière, CHU Nîmes, Nîmes, France

OPEN ACCESS

Edited by:

Maria Teresa Mascellino, Sapienza University of Rome, Italy

Reviewed by:

Sabrina John Moyo, University of Bergen, Norway Wonkeun Song, Hallym University, South Korea Guido Vincent Bloemberg, University of Zurich, Switzerland

*Correspondence:

Jean-Philippe Lavigne jean.philippe.lavigne@chu-nimes.fr

Specialty section:

This article was submitted to Clinical Microbiology, a section of the journal Frontiers in Cellular and Infection Microbiology

Received: 16 May 2020 Accepted: 06 August 2020 Published: 11 September 2020

Citation:

Durand C, Boudet A, Lavigne J-P and Pantel A (2020) Evaluation of Two Methods for the Detection of Third Generation Cephalosporins Resistant Enterobacterales Directly From Positive Blood Cultures. Front. Cell. Infect. Microbiol. 10:491. doi: 10.3389/fcimb.2020.00491 Due to the importance of a rapid determination of patients infected by multidrug resistant bacteria, we evaluated two rapid diagnostic tests for the detection of third-generation cephalosporins (3GC)-resistant Enterobacterales directly from positive blood cultures within 1 h: BL-REDTM (electrochemical method) and β -LACTATM test (chromogenic method). A panel of 150 clinical strains characterized for their resistance profiles (e.g., penicillinases, extended-spectrum beta-lactamases (ESBLs), overproduction of cephalosporinase, carbapenemases, impermeability) was tested. Approximately 100 CFU of each isolate was spiked into sterile blood culture bottles and incubated in a BD BACTECTM FX automated system (Becton Dickinson, USA). Positive blood cultures were examined to parallel testing using the BL-REDTM and β -LACTATM tests and conventional susceptibility method (disc diffusion following EUCAST recommendations). For all phenotypes combined, the sensitivity, specificity, positive predictive value, and negative predictive value in the detection of 3GC resistance were, respectively (i) with BL-REDTM: 45.7, 100, 100, and 54.2% and (ii) with β -LACTATM test: 52.2, 100, 100, and 56.9%. The positivity of tests allows to adapt antibiotic treatment whereas the negative result requires other tests. Moreover, these tests detect most Ambler class A-producing Enterobacterales (KPC, ESBL, extended-spectrum OXY) with sensitivities and specificities of 87.5 and 99% for BL-REDTM, respectively and both 100% for β-LACTATM test (47/47 isolates). These two rapid tests failed to detect AmpC overexpressed (sensitivities of 2.7% for BL-REDTM and 0% for β -LACTATM test) and Ambler class B-producing Enterobacterales (sensitivities of 40% for both tests) notably strains without ESBLs associated (sensitivities of 0% for both tests). BL-REDTM and β-LACTATM tests are easy-to-use and mainly attractive when a positive result is obtained notably to detect most of the Ambler class A-producing Enterobacterales in <1 h after the positivity of the blood culture, allowing a rapid adaptation of the antibiotic therapy in patients.

Keywords: bacteremia, BL-REDTM test, β -LACTATM test, blood cultures, extended-spectrum beta-lactamases, multidrug resistance, third-generation cephalosporins

INTRODUCTION

The spread of multidrug resistant (MDR) Enterobacterales is a major public health problem notably related to the misuse and overuse of antibiotics (Laxminarayan et al., 2013; World Health Organisation, 2014). Extended-spectrum cephalosporins and carbapenems are currently become the first line to treat community-acquired and nosocomial infections (Coque et al., 2008), involving a growing problem of the expansion of non-susceptibility of Enterobacterales to third-generation cephalosporins (3GC) and carbapenems. This increase is largely due to the importance of extended-spectrum β-lactamases (ESBL), but also to the overproduced chromosomal or plasmidmediated AmpC cephalosporinases and to the emergence of carbapenemases (Peirano and Pitout, 2019). The rapid optimization of antibiotic therapy, according to the organism and its resistance profile, is a major goal both for individual patients and for public health (Roca et al., 2015). Bloodstream infections (BSI) are 30% of admissions in ICU and is responsible for an increase in length of stay and a leading cause of mortality (close to 40%) (Kang et al., 2005). The importance of the early appropriate treatment is crucial in these pathologies, given the linear increase in the risk of mortality with each hour for which administration is delayed (Ibrahim et al., 2000; Ferrer et al., 2014). However, conventional microbiological methods are time consuming requiring 12-48 h to identify the causal microorganism and provide an antimicrobial susceptibility testing (AST). Over the last couple of years, the improvement of available rapid diagnostic tests directly from positive blood cultures has changed approaches for identification and AST (Maugeri et al., 2019). These tests have the advantage of providing fast results, optimizing antibiotic therapy, improving survival, and reducing the length of hospitalization (Perez et al., 2013). These new techniques promote the proper use of antibiotics by limiting the antimicrobial resistance and decreasing the medico-economic impact by (de)-escalation of the antibiotic therapy and reduction of hospital stay (Farfour et al., 2019). They include molecular biology tests identifying the MDR encoding genes (pathogen-specific real-time or multiplex PCR), spectrometry assays detecting *β*-lactamase hydrolytic activity (MALDI TOF-MS), biochemical tests, β-lactamases inhibitorbased tests, chromogenic tests and electrochemical test (Buchan et al., 2013; Renvoise et al., 2013; Bogaerts et al., 2016; Salimnia et al., 2016; Faron et al., 2017; Pantel et al., 2018).

The aim of the present study was to evaluate the performance of: (i) a new electrochemical test, BL-REDTM (Beta-Lactamase Rapid Electrochemical Detection, CORIS BioConcept, Belgium); and (ii) a colorimetric test, β -LACTATM test (Bio-Rad, Marnes la Coquette, France), for the detection of 3GC-resistance on a panel of Enterobacterales directly from positive blood cultures within 1 h.

MATERIALS AND METHODS

Bacterial Panel

A panel of 150 Enterobacterales isolates from our regional MDR Gram-negative Bacilli Reference Lab (CARB-LR group) in the

Occitanie region was tested (Pantel et al., 2014; Robert et al., 2014). The isolates were included with the following distribution: *Escherichia coli* (n = 62), *Klebsiella pneumoniae* (n = 29), *Enterobacter cloacae* (n = 20), *Klebsiella aerogenes* (n = 9), *Proteus mirabilis* (n = 7), *Morganella morganii* (n = 6), *Klebsiella oxytoca* (n = 5), *Citrobacter freundii* (n = 5), *Citrobacter koseri* (n = 2), *Providencia rettgeri* (n = 2), *Serratia marcescens* (n = 2), and *Hafnia alvei* (n = 1).

Different β -lactam resistance profiles were selected: (i) 3GCsusceptible isolates (n = 50), (ii) ESBL producers (n =41, including three which also overproduced AmpC), (iii) chromosomal-hyperproduced (n = 36) or plasmid-mediated (n =10) AmpC producers, (iv) extended-spectrum OXY (n =3), and (v) carbapenemase producers (n = 30, including 12 ESBL, 4 chromosomal AmpC, and 1 plasmid-mediated AmpC). Twenty-two isolates producing ESBL, high level AmpC or OXY were resistant to ertapenem by membrane permeability alteration. The presence of ESBLs, plasmid mediated AmpC, and carbapenemases genes was previously confirmed by specific PCR and sequence analysis (Perez-Perez and Hanson, 2002; Pitout et al., 2004, 2007; Poirel et al., 2011). All characteristics of the strains are summarized in **Table 1**.

Sample Preparation

For each isolate, ~ 100 CFU was spiked into sterile blood cultures [BD BACTECTM Plus Aerobic/F and BD BACTECTM Lytic/10 Anaerobic/F (BD Diagnostics, Le Pont de Claix, France)] containing 10 mL of fresh blood. Blood cultures were incubated in a BACTECTM FX automated blood culture device until they flagged positive for microbial growth. All positive blood cultures were divided into three samples: 1 mL for the conventional comparator methods, 1 mL for the β -LACTATM test (Bio-Rad, Marnes-la-Coquette, France) and 0.5 mL for the BL-REDTM test (Coris BioConcept, Gembloux, Belgium).

Conventional Comparator Method

The positive blood cultures were subcultured on blood agar (bioMérieux, Marcy l'Etoile, France) for 18–24 h at 37 \pm 2°C. Species-level identification was performed on isolated colonies on blood agar by the mass spectrometry (Vitek[®] MS, bioMérieux). AST was determined by the agar disc-diffusion method according to the EUCAST-SFM 2019 recommendations (www.eucast.org). The results of disc diffusion were used as the comparator for the 3GC resistance of each isolate (determined by the agar disc-diffusion size), using notably: ceftazidime 10 µg and cefotaxime 5 µg discs.

BL-RED[™] Test V1.0

All blood cultures were tested after the blood culture flagged positive as recommended by the manufacturer. For aerobic blood culture, a 40 μ L sample was mixed with 40 μ L of reagent and incubated 1 h at 37 \pm 2°C. For anaerobic blood culture, a pre-analysis was needed: a 500 μ L sample was centrifuged one time 2 min at 6,000 g; the pellet was re-suspended into 500 μ L of NaCl solution (0.9%) and re-centrifuged 2 min at 6,000 g; this second pellet was suspended into 40 μ L of buffer and 40 μ L of reagent; the mix was incubated 1 h at 37 \pm

TABLE 1 | Characteristics of the studied Enterobacterales isolates.

Resista	ance profile (no. of strains)	Species (no. of strains)	β -lactamase content (no. of strains)	
3GC susceptible (50)	No resistance to β -lactams (19)	<i>E. coli</i> (18)	None (18)	
		P. mirabilis (1)	None (1)	
	Penicillinases (22)	<i>E. coli</i> (17)	TEM-1 (10)	
			Inhibitor-resistant TEM (7)	
		K. pneumoniae (3)	SHV-1 (2)	
			Inhibitor-resistant TEM (1)	
		C. koseri (2)	CKO (2)	
	Cephalosporinases (9)	M. morganii (3)	Low-level AmpC (3)	
		P. rettgeri (2)	Low-level AmpC (2)	
		K. aerogenes (1)	Low-level AmpC (1)	
		E. cloacae (1)	Low-level AmpC (1)	
			Low-level AmpC (1)	
			Low-level AmpC (1)	
-lactamases conferring	ESBL (29)		CTX-M-group 1 (7)	
esistance to 3GC (70)	()		CTX-M-group 1 + CTX-M-group 8 (5)	
			ан ан ()	
		D mirabilia (E)	CTX-M-group 9 (4)	
		P. Miradilis (5)	CTX-M-group 1 (3)	
			CTX-M-group 9 (1)	
			CTX-M no group typed (1)	
		K. pneumoniae (4)	CTX-M-group $1 + $ SHV-1 (2)	
			SHV-5 (1)	
			CTX-M-group 1 + CTX-M-group 8 + SHV-1 (1)	
		E. cloacae (3)	CTX-M-group 1 + High-level AmpC (1)	
			CTX-M-group 9 + High-level AmpC (1)	
			TEM-24 + High-level AmpC (1)	
		K. oxytoca (1)	CTX-M-group 8 + OXY-1 (1)	
	Chromosomal overproduced cephalosporinases (29)	E. cloacae (9)	High-level AmpC (9)	
		E. coli (6)	High-level AmpC (6)	
		K. aerogenes (6)	High-level AmpC (6)	
		C. freundii (4)	High-level AmpC (4)	
		M. morganii (3)	High-level AmpC (3)	
		<i>H. alvei</i> (1)	High-level AmpC (1)	
	Plasmid cephalosporinases (9)	K. pneumoniae (9)	DHA-1 + SHV-1 (9)	
	Chromosomal overproduced penicillinases (3)	K. oxytoca (3)	High-level OXY-1 (3)	
arbapenemases (30)	Class A carbapenemase (6)	nce to β-lactams (19) E. coli (18) P. mirabilis (1) E. coli (17) K. pneumoniae (3) C. koseri (2) M. morganii (3) P. rettgeri (2) K. aerogenes (1) E. cloacae (1) E. cloacae (1) E. coli (1) S. marcescens (1) E. coli (16) P. mirabilis (5) K. pneumoniae (4) E. cloacae (3) K. oxytoca (1) E. cloacae (3) K. oxytoca (1) E. cloacae (9) E. coli (6) K. aerogenes (6) C. freundii (4) M. morganii (3) H. alvei (1) k. pneumoniae (5) probapenemase (10) E. cloacae (1) E. cloacae (1) E. cloacae (1) E. cloacae (1) E. cloacae (2) E. coli (1) K. pneumoniae (2) C. freundii (1) P. mirabilis (1)	KPC-2 + SHV-1 (4)	
			KPC-2 + CTX-M-group 1 + SHV-1 (1)	
		E. cloacae (1)	IMI-1 + Low-level AmpC (1)	
	Class B carbapenemase (10)	E. cloacae (4)	VIM-1 + High-level AmpC (3)	
			VIM-1 + CTX-M-group 9 (1)	
		E. coli (2)	NDM-1 + CTX-M-group 9 (1)	
			NDM-1 + DHA (1)	
		K. pneumoniae (2)	NDM-1 + CTX-M-group 1 + SHV-1 (2)	
			VIM-1 + High-level AmpC (1)	
			NDM-1 (1)	
	Class D carbapenemase (12)		OXA-48 + SHV-1 (2)	
		рноанноннае (4)	OXA-40 + SHV-1(2) OXA-48 + CTX-M-group 1 + SHV-1 (2)	
		K aprogramas (0)		
			OXA-48 + Low-level AmpC (2)	
		E. Cloacae (2)	OXA-48 + CTX-M-group 1 + Low-level AmpC	

(Continued)

TABLE 1 | Continued

Resistance profile (no. of strains)	Species (no. of strains)	β -lactamase content (no. of strains)
	E. coli (2)	OXA-48 (1)
		OXA-181 + CTX-M-group 1 (1)
	K. oxytoca (1)	OXA-48 + OXY-1 (1)
	S. marcescens (1)	OXA-48 + Low-level AmpC (1)
Class B + D carbapenemase (2)	K. pneumoniae (2)	NDM-1 + OXA-48 + CTX-M-group 1 + SHV-1 (2)

TABLE 2 | Analytical performance of assays for detection of 3GC-resistant Enterobacterales isolates in positive blood cultures.

Method	Result	3GC cephalosporins		Sensitivity (95% CI) Specificity (95% CI) PPV (95% CI)			NPV (95% CI)	Accuracy (95% CI) Youden	
		Resistant (n = 92)	Susceptible (n = 58)	_					
BL-RED [™]	Positive	43	0	46.7% (36.9–56.9)	100% (93.8–100)	100% (91.8–10)	54.2% (44.8–63.3)	67.3%	46.7%
	Negative	49	58						
β -LACTA TM	Positive	48	0	52.2% (42.1–62.1)	100% (93.8–100)	100% (92.6–100)	56.9% (47.2–66.1)	70.7%	52.2%
	Negative	44	58						

PPV, positive predictive value; NPV, negative predictive value.

 $2^\circ C.$ After this incubation, 20 μL of the mix was placed on the ceramic electrode and after in the sensor. The reagent includes a "false" 3GC substrate which is hydrolysed when there is a 3GC β -lactamase, this frees an electro-conductive product. The intensity of the electrochemical signal [measured by DropSens (Metrohm, Villebon Courtaboeuf, France)] is posted in nanoAmpere (nA). The threshold of positivity is 80 nA for blood culture.

β -LACTATM Test

The protocol B of the manufacturer recommendations was followed for the treatment of the positive blood cultures. Three different centrifugations and different solutions (Triton 0.01%, NaCl solution) were used to obtain the final bacterial pellet. After using the reagents provided in the kit on the bacterial pellet, we can read visually the hydrolysis of the chromogenic substrate (HMRZ-86, cephalosporin) from yellow to red in 15 min for a positive test. No change in color was considered a negative result (no hydrolysis of HMRZ-86).

Statistical Analysis

The sensitivity, specificity, negative and positive predictive value (NPV and PPV, respectively) were calculated and comparison with molecular characterization which served as the reference standard. 3GC-sensitive isolate with a negative test and 3GC-resistant isolate with a positive test are considered correct. Conversely, 3GC-sensitive isolate with a negative test are considered incorrect. Additionally, the 95% confidence intervals (CIs), accuracy and the Youden index were calculated.

RESULTS

Evaluation on 3GC-Susceptible Isolates

Among our panel, 50 isolates had no 3GC resistance. These isolates harbored either a penicillinase (CKO, TEM-1, SHV-1, inhibitor-resistant TEM) or a naturally chromosomally-encoded inducible AmpC cephalosporinase. The two assays confirmed the absence of 3GC-resistance in all strains: 0 nA for most of the strains with BL-REDTM, yellow color with β -LACTATM test.

Of note, two isolates (one wild-type *E. coli* and one TEM-1-producing *E. coli*) presented a very low intensity using BL-REDTM at 10 nA. However, they were still considered negative because the positive threshold for blood culture was established at 80 nA.

Global Detection of 3GC-Resistance

Concerning all isolates belonging to the tested panel, the performance of the BL-REDTM and β -LACTATM tests were, respectively: sensitivity, 46.7% (CI 36.9–56.9%) and 52.2% (CI 42.1–62.1%); specificity, both 100% (CI 93.8–100%); NPV, 54.2% (CI 44.8–63.3%), and 56.9% (CI 47.2–66.1%) and PPV, both 100% (CI 91.8–100%; 92.6–100%). The Youden index was 46.7 and 52.2%, respectively (**Table 2**).

Evaluation on Ambler Class a β -Lactamase Producers

When considering only the Ambler class A ESBL (without carbapenemase associated), 79.3% (23/29) of the isolates were positive using the BL-REDTM, with a signal >80 nA [100–1,220] (**Table S1**). Six isolates were negative: two CTX-M-group 1-producing *P. mirabilis*, one CTX-M-group 9-producing *E. coli*, one CTX-M-group 9 and 8-producing *E. coli*, one CTX-M-group 9-producing *E. cloacae* with altered membrane permeability, and one SHV-5-producing *K. pneumoniae* with altered membrane

Method	Result -	Class A b-lactamase*		Sensitivity (95% CI) Specificity (95% CI) PPV (95% CI)			NPV (95% CI)	Accuracy (95% CI) Youden	
			Susceptible (n = 113)						
BL-RED [™]	Positive	31	1	83.8% (68.9–92.4)	99.1% (95.2–99.8)	96.9% (84.3–99.5)	94.9% (89.4–97.7)	95.3%	82.9%
	Negative	6	112						
$\beta\text{-LACTA}^{\text{TM}}$	Positive Negative	37 0	0 113	100% (90.6–100)	100% (96.7–100)	100% (90.6–100)	100% (96.7–100)	100%	100%

TABLE 3 Analytical performance of assays for detection of Ambler Class A β-lactamases-producing Enterobacterales isolates in positive blood cultures.

*IMI-producing E. cloacae strain susceptible to 3GC was excluded.

PPV, Positive Predictive Value; NPV, Negative Predictive Value.

permeability. All these isolates (n = 29) were positive using the β -LACTATM test that shows a sensitivity of ESBL detection of 100%. Interestingly, the two techniques detected the OXY-hyperproducing *K. oxytoca* isolates with low outer membrane permeability (n = 3).

Concerning Ambler class A carbapenemases, both tests were able to detect the KPC-producing isolates (n = 5) whatever the presence or absence of an ESBL (**Table S1**). Negative results were obtained with IMI-1-producing *E. cloacae* strain that is resistant to carbapenems but remains susceptible to 3GC.

The analysis of the different previous results suggested that the two methods evaluated in this study had an interest in the detection of Ambler class A producers (not restricted to ESBL producers). When considering all Ambler class A producers (excluding IMI producers), the performance of BL-REDTM and β -LACTATM were, respectively: sensitivity, 87.5% (CI 75.3–94.1%) and 100% (CI 92.6–100%); specificity, 99.0% (CI 94.7–99.8%) and 100% (CI 96.4–100%); PPV, 97.7% (CI 87.9–99.6%) and 100% (CI 96.4–100%) and NPV, 94.4% (CI 88.3–97.4%) and 100% (CI 96.4–100%) (**Table 3**).

Evaluation on Ambler Class C β-Lactamase Producers

All AmpC-overproducing isolates (n = 29, excluding ESBL co-producers) tested were negative for the two techniques (**Table S2**). The alteration of outer membrane permeability (seven strains) did not influence the results. One isolate (*C. freundii*) presented a high intensity using BL-REDTM at 60 nA but under the positive threshold. The antibiogram analysis performed after blood cultures confirmed the resistance to 3GC and the presence of the hyperproduction of AmpC in all isolates.

Among plasmid-mediated AmpC producers (n = 10), only one strain (*K. pneumoniae* DHA-1 with decreased membrane permeability) was positive with BL-REDTM test (480 nA).

Evaluation on Ambler Class B and D β-Lactamase Producers


Despite the high level of 3GC resistance, Ambler class B carbapenemase producers were only detected (6/12) in both tests when the isolates were also ESBL co-producers.

Similarly, OXA-48 carbapenemases were only detected (7/14) by the two tests when the isolates were ESBL co-producers.

DISCUSSION

Enterobacterales are the most important etiologies of community-onset and hospital-acquired BSI (Laupland and Church, 2014). The diffusion of 3GC-resistance among Enterobacterales, resulting from the spread of ESBL- and carbapenemase-producing isolates represents a serious problem worldwide. This is particularly true for life-threatening infections such as BSI, for which inappropriate first-line antibiotic therapy has a dramatic impact on short-term mortality (Kumar et al., 2009; Zilberberg et al., 2014).

Rapid and reliable identification of MDR Enterobacterales in the clinical microbiology laboratory is essential to effective infection control. In this study, we evaluated two rapid and easy-to-use methods: BL-RED^{TM} and β -LACTA^{TM} tests. If the β-LACTATM test has been used for few years on bacterial strains and directly from blood cultures, the present study is the first to evaluate the recently commercialized BL-REDTM test v1.0. Both tests allowed early detection of most Ambler class A β-lactamases conferring 3GC resistance (e.g., ESBL, OXY, and KPC) to Enterobacterales (Table 3) in <1 h after positivity of the blood culture with interesting performance (sensitivity and specificity to 83.8 and 99.1% for BL-REDTM test and both 100% for β -LACTATM test) suggesting that the two tests may provide useful therapeutic guidance. However, when testing challenged 3GC-resistant isolates non-ESBL producers, the two tests showed poor performance. Of note, a very low detection rate of plasmid-mediated AmpC or AmpCoverproducing isolates could be observed: 2.4% for BL-REDTM and 0% for β -LACTATM test. This limitation has been previously published concerning β -LACTATM test, when testing with subculture on agar plates (Renvoise et al., 2013; Morosini et al., 2014; Compain et al., 2015). The lack of sensitivity for the detection of high level AmpC must be balanced against the fact that ESBL producers are predominant among 3GC resistant Enterobacterales in European countries, notably among BSI (European Antimicrobial Resistance Surveillance Network (EARS-Net), 2018). Moreover, early adaptation by cefepime could be proposed when naturally AmpC producer was identified by MALDI-TOF MS, as previously suggested (Mizrahi et al., 2018). Another limitation of both tests was the failure to identify metallo-β-lactamase-producing isolates directly in blood cultures, whereas β-LACTATM test carried

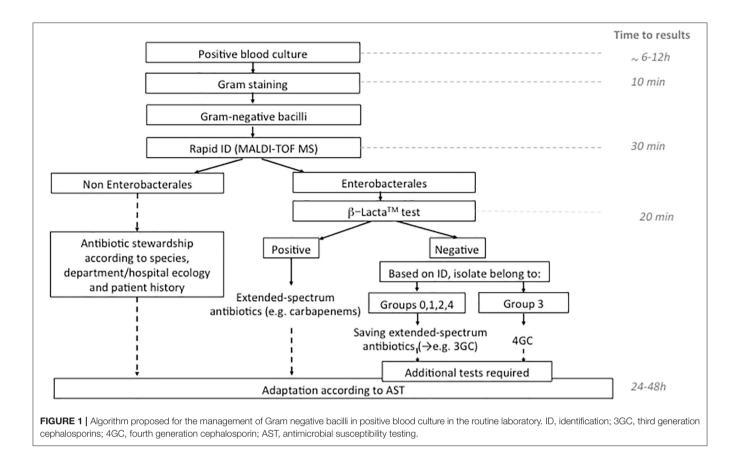

out with VIM-1-producers on agar plates showed acceptable sensitivity (Morosini et al., 2014). The lack of detection of metallo- β -lactamase-producing Enterobacterales isolates have been previously reported with phenotypical tests applied on blood culture, the zinc concentration in the medium being essential for an efficient detection of Ambler class B enzymes (Dortet et al., 2014; Pantel et al., 2018). Concerning class D carbapenemase producers, the isolates are frequently associated with ESBL positive producers (Bakthavatchalam et al., 2016) allowing their detection with both evaluated tests (7/7, 100%). Without the presence of ESBL-carrying plasmids, the tests failed to identify this resistance mechanism. To date, molecular methods remain the faster tool with a higher sensitivity to detect the carbapenemase-producing Enterobacterales (Rood and Li, 2017).

If no study has previously evaluated the BL-REDTM test, some previous works were performed on the use of β -LACTATM test combined with MALDI-TOF MS from positive blood cultures (Compain et al., 2015; Walewski et al., 2015; Hasso et al., 2017; Mizrahi et al., 2018). The first study conducted in 2015 by Compain et al. evaluated the β -LACTATM test on 3-h subcultures. 84.8% of the 33 blood culture isolates resistant to 3GC were correctly detected Compain et al. (2015). Walewski et al. (2015) used the test on bacterial pellets from the blood culture broths after treatment by saponin followed by two washes, with 95.7% of sensitivity and 100% of specificity for identifying ESBLproducing Enterobacterales. Hasso et al. (2017) evaluated the accuracy of the β -LACTATM test for rapid detection of ESBLproducing *E. coli* and *Klebsiella* spp. from smudge plates prepared with positive blood cultures. The authors noted a sensitivity and specificity of 100 and 97.8%, respectively. Finally, the last two studies carried out in France by Mizrahi et al. (2018) and Farfour et al. (2019) obtained similar results: excellent performance regarding ESBL-producers and lack of sensitivity for Group 3 Enterobacterales. All results underlined the interest of these rapid techniques to be used effectively in daily laboratory practice.

The use of rapid phenotypic susceptibility techniques have demonstrated their interest in an optimized antimicrobial treatment (Schneider et al., 2019). Globally, our study shows excellent PPV (both 100%) but poor NPV (46.7 and 52.2%) for the BL-REDTM and β -LACTATM tests, respectively. This reinforces the idea to only manage positive results with these assays and adapt antibiotic treatment whereas the negative results need to use additional tests. In this way, we suggest an algorithm incorporating MALDI-TOF MS combined with β-LACTATM test for optimized Ambler class A producers detection with the advantage of being rapid, affordable, and simple (Figure 1). This algorithm could be accurate to give in <2 h a decisive orientation in the antibiotic management of sepsis or septic shock patients, particularly when the β -LACTATM test is positive. An early adaptation by 4GC could be proposed when naturally AmpC producer was identified by MALDI-TOF MS with a negative β-LACTATM test. However, its use must be debated. Recently, Dépret et al. (2018) evaluated the rapid identification

using MALDI-TOF MS results only vs. the combined MALDI-TOF MS and the β -LACTATM test in a method close to our algorithm. No impact on the current standard choice of therapy in terms of escalation/de-escalation or in reduction of carbapenems prescription was observed. Indeed, only 9% of the positive blood cultures in this study were confirmed as ESBL-producing, decreasing the interest of the techniques. This suggests that the use of a rapid test for Ambler class A-producers detection would be useful in a population with high prevalence of these enzymes and/or when treatment choices are not made by infectious disease specialists. In contrast, Mizrahi et al. (2018) have demonstrated the crucial role played by rapid detection of 3GC resistance with β -LACTATM. In this study, 75% (21/28) of the patients with BSI involving 3GC-resistant bacteria received a non-adapted first-line treatment. β -LACTATM test performed on blood culture significantly reduced the delay for treatment adaptation (28.1 h) and patient isolation (35 h). In 2019, Farfour et al. (2019) observed the same impact of rapid diagnostic tests on the management of BSI. Patients receiving a rapid strategy (rapid identification and 3GC-susceptibility testing) received more frequently an effective and appropriate antibiotic therapy than patients receiving conventional strategy on the first day of BSI diagnosis. Studies in settings where the prevalence of Ambler class A producers is high must be done to definitively conclude on the interest of the rapid technique.

In conclusion, the two evaluated assays (BL-REDTM and β -LACTATM tests) showed excellent PPV, at least for the detection of the common Ambler class A family conferring 3GC-resistance and are suitable for the routine microbiology laboratory, decreasing the time to detection, and the need to refer these isolates to molecular biology confirmation. When a positive result is detected, these techniques permitted a quick and easy detection of the presence of ESBL-producing Enterobacterales in a clinical sample. They may help clinicians to guide appropriate antimicrobial therapy in septic patients and presumably improve their prognosis in a country or region where the Ambler class A family is particularly prevalent.

REFERENCES

- Bakthavatchalam, Y. D., Anandan, S., and Veeraraghavan, B. (2016). Laboratory detection and clinical implication of oxacillinase-48 like carbapenamse: the hidden threat. J. Glob. Infect. Dis. 8, 41–50. doi: 10.4103/0974-777X.176149
- Bogaerts, P., Yunus, S., Massart, M., Huang, T. D., and Glupczynski, Y. (2016). Evaluation of the BYG carba test, a new electrochemical assay for rapid laboratory detection of carbapenemase-producing *Enterobacteriaceae. J. Clin. Microbiol.* 54, 349–358. doi: 10.1128/JCM.02404-15
- Buchan, B. W., Ginocchio, C. C., Manii, R., Cavagnolo, R., Pancholi, P., Swyers, L., et al. (2013). Multiplex identification of Gram-positive bacteria and resistance determinants directly from positive blood culture broths: evaluation of an automated microarray-based nucleic acid test. *PLoS Med.* 10:e1001478. doi: 10.1371/journal.pmed.1001478
- Compain, F., Bensekhri, H., Rostane, H., Mainardi, J. L., and Lavollay, M. (2015). β LACTA test for rapid detection of *Enterobacteriaceae* resistant to third-generation cephalosporins from positive blood cultures using briefly incubated solid medium cultures. *J. Med. Microbiol.* 64, 1256–1259. doi: 10.1099/jmm.0.000157

DATA AVAILABILITY STATEMENT

All datasets presented in this study are included in the article/**Supplementary Material**.

AUTHOR CONTRIBUTIONS

CD has performed experiments, statistical analysis, and wrote the manuscript. AB has performed experiments and critically reviewed the manuscript. J-PL and AP have conceived the study, discussed the results, and wrote the manuscript. All authors contributed to the article and approved the submitted version.

FUNDING

This work was supported by the University Hospital of Nîmes.

ACKNOWLEDGMENTS

We thank Coris BioConcept, Gembloux, Belgium, for providing BL-REDTM v.1.0 kit with test reagents. This company had no role in the study design, data collection or interpretation of the results.

AP and J-PL belong to the FHU INCh (Federation Hospitalo Universitaire Infections Chroniques, Aviesan). We thank the Nîmes University Hospital for its structural, human and financial support through the award obtained by our team during the internal call for tenders \ll Thématiques phares \gg .

We thank Sarah Kabani for her editing assistance.

This work has been selected at the 30th European Congress of Clinical Microbiology and Infectious Diseases, Paris (France) April, 2020 and presented in the Abstract book (8551).

SUPPLEMENTARY MATERIAL

The Supplementary Material for this article can be found online at: https://www.frontiersin.org/articles/10.3389/fcimb. 2020.00491/full#supplementary-material

- Coque, T. M., Baquero, F., and Canton, R. (2008). Increasing prevalence of ESBL-producing *Enterobacteriaceae* in Europe. *Euro Surveill*. 13:19044. doi: 10.2807/ese.13.47.19044-en
- Dépret, F., Aubry, A., Fournier, A., Charles-Nelson, A., Katsahian, S., Compain, F., et al. (2018). β LACTA testing may not improve treatment decisions made with MALDI-TOF MS-informed antimicrobial stewardship advice for patients with Gram-negative bacteraemia: a prospective comparative study. *J. Med. Microbiol.* 67, 183–189. doi: 10.1099/jmm.0.000665
- Dortet, L., Bréchard, L., Poirel, L., and Nordmann, P. (2014). Impact of the isolation medium for detection of carbapenemase-producing *Enterobacteriaceae* using an updated version of the Carba NP test. J. Med. Microbiol. 63, 772–776. doi: 10.1099/jmm.0.071340-0
- European Antimicrobial Resistance Surveillance Network (EARS-Net) (2018). Annual Report of the European Antimicrobial Resistance Surveillance Network (EARS-Net). Available online at: https://www.ecdc.europa.eu/ (accessed May 10, 2020).
- Farfour, E., Si Larbi, A. G., Cardot, E., Limousin, L., Mathonnet, D., Cahen, P., et al. (2019). Impact of rapid diagnostic tests on the management of patients presenting with *Enterobacteriaceae* bacteremia. *Med. Mal. Infect.* 49, 202–207. doi: 10.1016/j.medmal.2018.11.015

- Faron, M. L., Buchan, B. W., and Ledeboer, N. A. (2017). Matrix-Assisted Laser Desorption Ionization – Time of Flight Mass Spectrometry for use with positive blood cultures: methodology, performance and optimization. *J. Clin. Microbiol.* 55, 3328–3338. doi: 10.1128/JCM.00868-17
- Ferrer, R., Martin-Loeches, I., Phillips, G., Osborn, T. M., Townsend, S., Dellinger, R. P., et al. (2014). Empiric antibiotic treatment reduces mortality in severe sepsis and septic shock from the first hour: results from a guidelinebased performance improvement program. *Crit. Care Med.* 42, 1749–1755. doi: 10.1097/CCM.00000000000330
- Hasso, M., Porter, V., and Simor, A. E. (2017). Evaluation of the β-Lacta Test for detection of Extended-Spectrum-β-Lactamase (ESBL)-producing organisms directly from positive blood cultures by use of smudge plates. J. Clin. Microbiol. 55, 3560–3562. doi: 10.1128/JCM.01354-17
- Ibrahim, E. H., Sherman, G., Ward, S., Fraser, V. J., and Kollef, M. H. (2000). The influence of inadequate antimicrobial treatment of bloodstream infections on patient outcomes in the ICU setting. *Chest* 118, 146–155. doi: 10.1378/chest.118.1.146
- Kang, C., Kim, S., Park, W., Lee, K., Kim, H., Kim, E., et al. (2005). Bloodstream infections caused by antibiotic-resistant Gram-negative bacilli : risk factors for mortality and impact of inappropriate initial antimicrobial therapy on outcome. *Antimicrob. Agents Chemother.* 49, 760–766. doi: 10.1128/AAC.49.2.760-766.2005
- Kumar, A., Ellis, P., Arabi, Y., Roberts, D., Light, B., Parrillo, J. E., et al. (2009). Initiation of inappropriate antimicrobial therapy results in a fivefold reduction of survival in human septic shock. *Chest* 136,1237–1248. doi: 10.1378/chest.09-0087
- Laupland, K. B., and Church, D. L. (2014). Population-based epidemiology and microbiology of community-onset bloodstream infections. *Clin. Microbiol. Rev.* 27, 647–664. doi: 10.1128/CMR.00002-14
- Laxminarayan, R., Duse, A., Wattal, C., Zaidi, A. K. M., Wertheim, H. F. L., Sumpradit, N., et al. (2013). Antibiotic resistance-the need for global solutions. *Lancet Infect. Dis.* 13, 1057–1098. doi: 10.1016/S1473-3099(13)70318-9
- Maugeri, G., Lychko, I., Sobral, R., and Roque, A. C. A. (2019). Identification and antibiotic-susceptibility profiling of infectious bacterial agents: a review of current and future trends. *Biotechnol. J.* 14:e1700750. doi: 10.1002/biot.201700750
- Mizrahi, A., Amzalag, J., Couzigou, C., Péan De Ponfilly, G., Pilmis, B., and Le Monnier, A. (2018). Clinical impact of rapid bacterial identification by MALDI-TOF MS combined with the bêta-LACTATM test on early antibiotic adaptation by an antimicrobial stewardship team in bloodstream infections. *Infect. Dis.* 50, 668–677. doi: 10.1080/23744235.2018.1458147
- Morosini, M. I., García-Castillo, M., Tato, M., Gijón, D., Valverde, A., Ruiz-Garbajosa, P., et al. (2014). Rapid detection of β -lactamase-hydrolyzing extended-spectrum cephalosporins in *Enterobacteriaceae* by use of the new chromogenic β Lacta test. *J. Clin. Microbiol.* 52, 1741–1744. doi: 10.1128/JCM.03614-13
- Pantel, A., Boutet-Dubois, A., Jean-Pierre, H., Marchandin, H., Sotto, A., and Lavigne, J. P. (2014). French regional surveillance program of carbapenemaseproducing Gram-negative bacilli: results from a 2-year period. *Eur. J. Clin. Microbiol. Infect. Dis.* 33, 2285–2292. doi: 10.1007/s10096-014-2189-5
- Pantel, A., Monier, J., and Lavigne, J. P. (2018). Performance of the Accelerate PhenoTM system for identification and antimicrobial susceptibility testing of a panel of multidrug-resistant Gram-negative bacilli directly from positive blood cultures. J. Antimicrob. Chemother. 73, 1546–1552. doi: 10.1093/jac/dky032
- Peirano, G., and Pitout, J. D. D. (2019). Extended-spectrum beta-lactamaseproducing *Enterobacteriaceae*: Update on molecular epidemiology and treatment options. *Drugs*. 79, 1529–1541. doi: 10.1007/s40265-019-01180-3
- Perez, K. K., Olsen, R. J., Musick, W. L., Cernoch, P. L., Davis, J. R., Land, G. A., et al. (2013). Integrating rapid pathogen identification and antimicrobial stewardship significantly decreases hospital costs. *Arch. Pathol. Lab. Med.* 137, 1247–1254. doi: 10.5858/arpa.2012-0651-OA

- Perez-Perez, F. J., and Hanson, N. D. (2002). Detection of plasmid-mediated AmpC beta-lactamase genes in clinical isolates by using multiplex PCR. J. Clin. Microbiol. 40, 2153–2162. doi: 10.1128/JCM.40.6.2153-2162.2002
- Pitout, J. D. D., Hamilton, N., Church, D. L., Nordmann, P., and Poirel, L. (2007). Development and clinical validation of a molecular diagnostic assay to detect CTX-M-type β -lactamases in *Enterobacteriaceae*. *Clin. Microbiol. Infect.* 13, 291–297. doi: 10.1111/j.1469-0691.2006.01645.x
- Pitout, J. D. D., Hossain, A., and Hanson, N. D. (2004). Phenotypic and molecular detection of CTX-M-β-lactamases produced by *Escherichia coli* and *Klebsiella* spp. J. Clin. Microbiol. 42, 5715–5721. doi: 10.1128/JCM.42.12.5715-5721.2004
- Poirel, L., Walsh, T. R., Cuvillier, V., and Nordmann, P. (2011). Multiplex PCR for detection of acquired carbapenemase genes. *Diagn. Microbiol. Infect. Dis.* 70, 119–123. doi: 10.1016/j.diagmicrobio.2010.12.002
- Renvoise, A., Decre, D., Amarsy-Guerle, R., Huang, T. D., Jost, C., Podglajen, I., et al. (2013). Evaluation of the βLacta test, a rapid test detecting resistance to third-generation cephalosporins in clinical strains of *Enterobacteriaceae. J. Clin. Microbiol.* 51,4012–4017. doi: 10.1128/JCM.01936-13
- Robert, J., Pantel, A., Mérens, A., Lavigne, J.-P., Nicolas-Chanoine, M.-H., Brieu, N., et al. (2014). Incidence rates of carbapenemase-producing *Enterobacteriaceae* clinical isolates in France: a prospective nationwide study in 2011-12. *J. Antimicrob. Chemother*. 69, 2706–2712. doi: 10.1093/jac/dku208
- Roca, I., Akova, M., Baquero, F., Carlet, J., Cavaleri, M., Coenen, S., et al. (2015). The global threat of antimicrobial resistance: science for intervention. *New Microbes New Infect.* 6, 22–29. doi: 10.1016/j.nmni.2015.02.007
- Rood, I. G. H., and Li, Q. (2017). Molecular detection of extended spectrum-β-lactamase- and carbapenemase-producing *Enterobacteriaceae* in a clinical setting. *Diagn. Microbiol. Infect. Dis.* 89, 245–250. doi:10.1016/j.diagmicrobio.2017.07.013
- Salimnia, H., Fairfax, M. R., Lephart, P. R., Schreckenberger, P., DesJarlais, S. M., Johnson, J. K., et al. (2016). Evaluation of the FilmArray blood culture identification panel: results of a multicenter controlled trial. *J. Clin. Microbiol.* 54, 687–698. doi: 10.1128/JCM.01679-15
- Schneider, J. G., Wood, J. B., Schmitt, B. H., Emery, C. L., Davis, T. E., Smith, N. W., et al. (2019). Susceptibility provision enhances effective de-escalation (SPEED): utilizing rapid phenotypic susceptibility testing in Gram-negative bloodstream infections and its potential clinical impact. J. Antimicrob. Chemother. 74, 116–123. doi: 10.1093/jac/dky531
- Walewski, V., Podglajen, I., Lefeuvre, P., Dutasta, F., Neuschwander, A., Tilouche, L., et al. (2015). Early detection with the β-LACTATM test of extended-spectrum β-lactamase-producing *Enterobacteriaceae* in blood cultures. *Diagn. Microbiol. Infect. Dis.* 83, 216–218. doi: 10.1016/j.diagmicrobio.2015.07.005
- World Health Organisation (2014). Antimicrobial Resistance: Global Report on Surveillance. Available online at: https://www.who.int/antimicrobialresistance/publications/surveillancereport/en/ (accessed May 10, 2020).
- Zilberberg, M. D., Shorr, A. F., Micek, S. T., Vazquez-Guillamet, C., and Kollef, M. H. (2014). Multi-drug resistance, inappropriate initial antibiotic therapy and mortality in Gram-negative severe sepsis and septic shock: a retrospective cohort study. *Crit. Care* 18, 1–13. doi: 10.1186/s13054-014-0596-8

Conflict of Interest: The authors declare that the research was conducted in the absence of any commercial or financial relationships that could be construed as a potential conflict of interest.

Copyright © 2020 Durand, Boudet, Lavigne and Pantel. This is an open-access article distributed under the terms of the Creative Commons Attribution License (CC BY). The use, distribution or reproduction in other forums is permitted, provided the original author(s) and the copyright owner(s) are credited and that the original publication in this journal is cited, in accordance with accepted academic practice. No use, distribution or reproduction is permitted which does not comply with these terms.